

1 Prepare the Reservoir:




Remove the reservoir from base and fill with lukewarm water. Replace the reservoir onto the base, pressing down firmly to ensure full seating.



2 Insert Tip:


Select the appropriate tip and insert into the unit. Press firmly until the tip clicks into place.

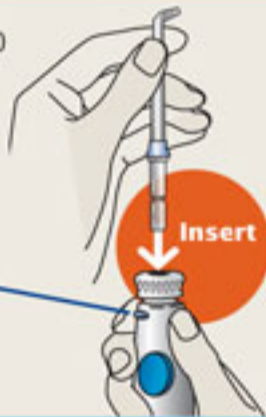


Caution: Be careful not to accidentally press the tip eject button during use.

3 Adjust Water Pressure:

For first time use, adjust the pressure dial on the base to pressure **Setting 1**. Gradually increase pressure as you prefer.

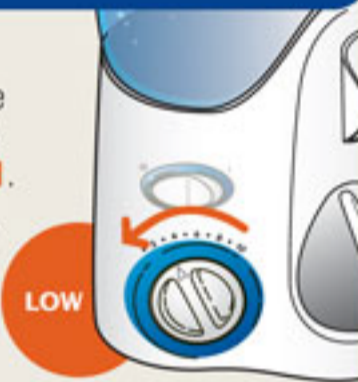

4 Handle Controls:

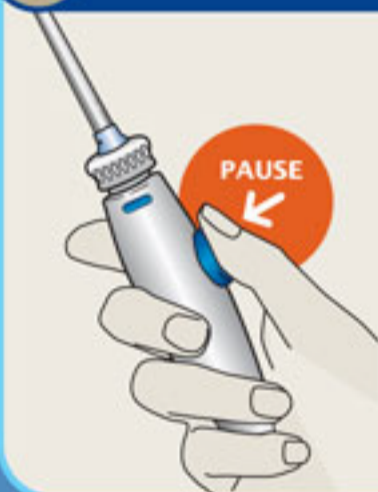

The **Waterpik Ultra** is equipped with a convenient pause button on the handle. You can temporarily stop water flow by pressing the pause button on the handle.



Direct water stream by turning the knob at the top of the handle.


5 ON/Off Procedure:



When tip is in your mouth, turn unit **ON** with your free hand. Close lips to prevent splashing.

Always **PAUSE** or turn unit **OFF** before removing tip from mouth.

6 Technique Tips for Best Results:

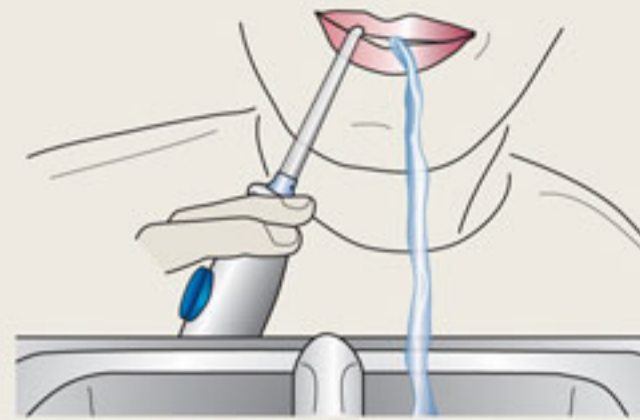

Place tip in mouth. Lean over sink and close lips enough to prevent splashing, while still allowing water to flow from mouth into the sink. Turn unit **ON**.




7 Technique Tips:

Aim the tip at the teeth at a 90 degree angle.



Starting with the back teeth, glide the tip along the gum line and pause briefly between teeth. Continue until all areas around and between teeth have been cleaned.

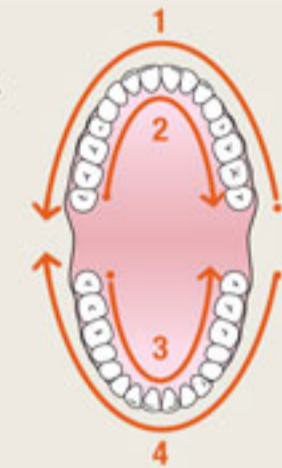

8 When Finished:

Turn the unit **OFF** and empty any liquid left in the reservoir.



Changing Tips:



To remove a tip, press the tip eject button on the side of the handle and pull tip out.

Caution: Be careful not to accidentally press the tip eject button during use.

NuDento International Limited

19/F., Yat Chau Building
262 Des Voeux Road Central
Hong Kong

Tel : (852) 3176 0561
Fax : (852) 3972 7278

www.waterpik.com.hk